


Sample Extractor

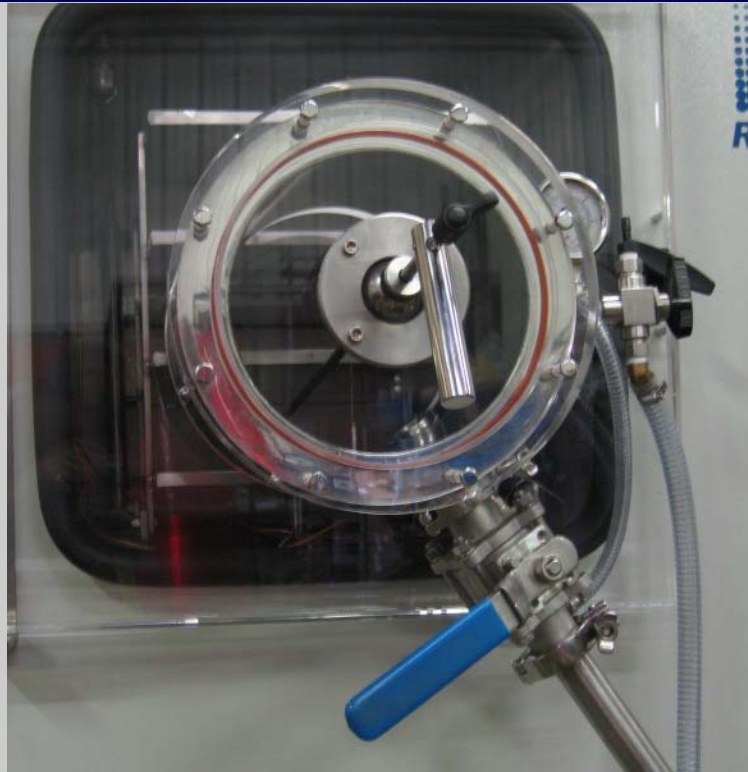

Standard Features	
Access Area	8" diameter
Reach	15" into chamber
Vial Sizes	Up to 20mm neck ID
Vial Accumulator	12", up to 5 each 10 ml vials
Backfill	Included
Stoppering	Included

Key Features:

- Full visibility
- Easy to use gripper
- Vacuum pump included
- 5 vial accumulator
- Backfill
- Stoppering



Millrock offers a unique solution to sample extraction. A 8" diameter opening offers unimpeded access to multiple shelves. The gripping mechanism reaches 15" into the chamber allowing vials to be easily selected from the sides and middle of the shelf.

The 12" vial accumulator speeds the collection process by providing room for one to five 10 ml vials to be selected for backfill and stoppering in each batch.

**MILLROCK TECHNOLOGY, Inc**

Kingston, NY 12401

845-339-5700

sales@millrocktech.com

www.millrocktech.com

